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
ROLLSTÄNDER ULTRA-LIGHT DUTY

ASSEMBLY INSTRUCTION

CART ULTRA-LIGHT DUTY

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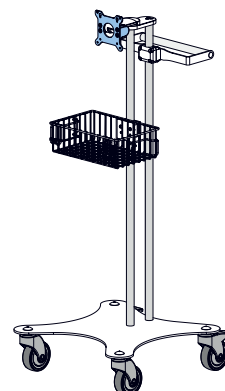Name: CIM Cart ultra-light duty
Part no. code of the cart: R-003
Part no. code of the monitor adaptation: A, FK, G, PH

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WARNING

Please check the weight of the ME equipment to be mounted. The permissible total weight is mentioned on the corresponding data sheet.

Failure to follow these warnings may result in damage to equipment or injury to people.

CIM med GmbH cannot be held responsible for any damage due to poor installation.


Product modifications, disassembling and reassembling and/or mounting any of the components in any other way than described below are permitted only after written approval by the CIM med GmbH construction and engineering department. CIM med GmbH cannot be held responsible for any damage or injury caused by product modifications and/ or amendments without prior permission from CIM med GmbH.

The additional instructions for use and the inspection plan are intended for end users.

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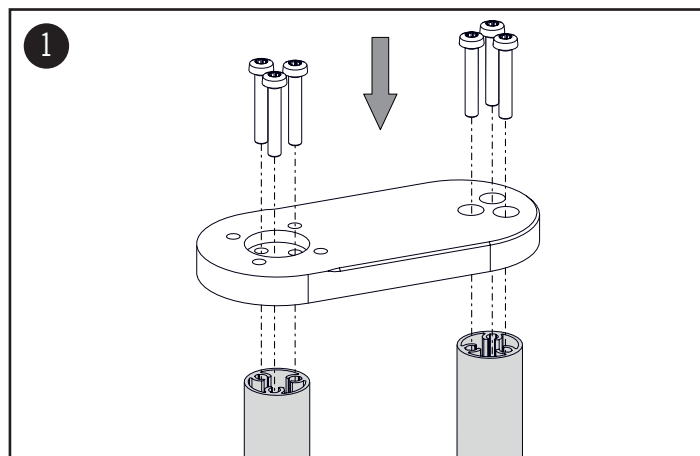
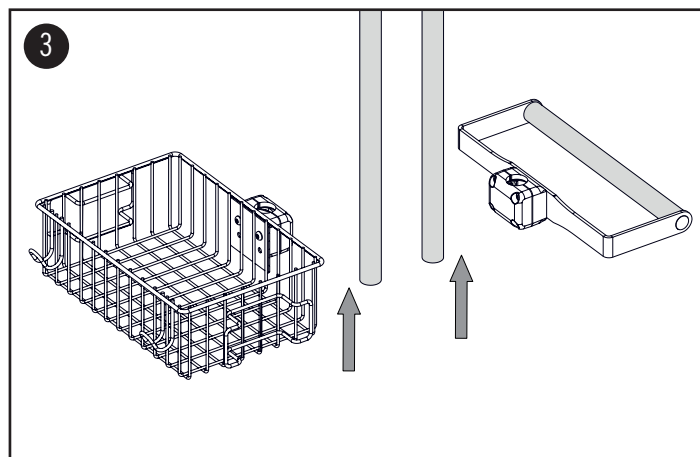
ASSEMBLY STEPS

ASSEMBLY STEPS


1. Mount the connection plate with the (6) M5x30 pan-head screws to the two poles.
2. Mount the monitor adaption with the (4) M5x20 countersunk screws.

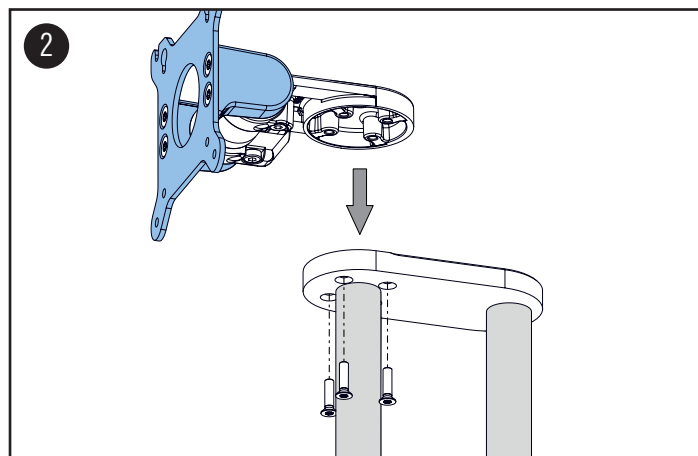3. Attach the accessory to the relevant pole at the desired height.

NOTE: Please refer to the corresponding assembly instruction of the accessory.

4. Mount the two poles to the base with (6) M5x30 pan-head screws. Tighten the screws securely.

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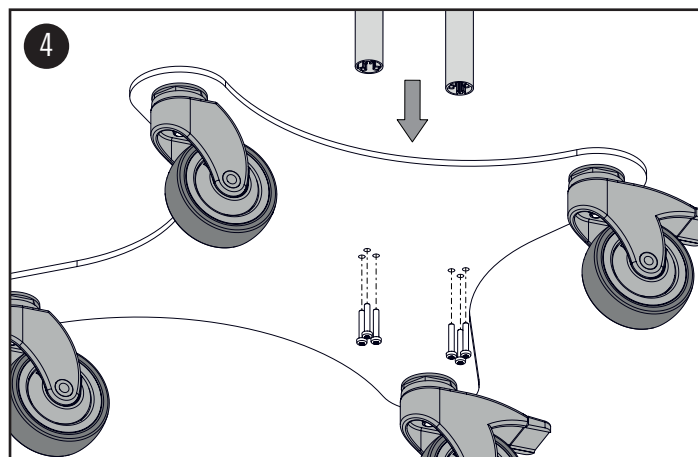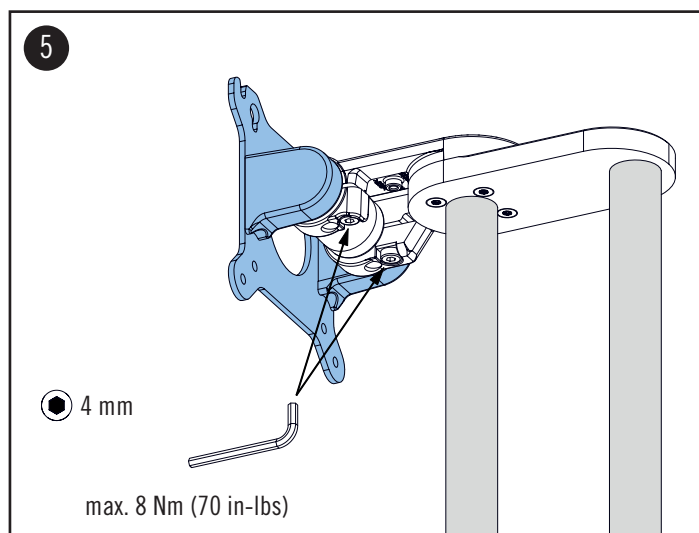
ASSEMBLY STEPS

TILT ADJUSTMENT

- To adjust the overall tilt tension, tighten or loosen the both hex screws. By tightening both of the screws the tilt function can be locked completely. Maximum permissible tightening torque is 8 Nm (70 in-lbs).

⚠ WARNING: In case of overtightening the screw the casting repectively the threads may break which leads to loss of function.

⚠ WARNING: When loosening both hex screws the ME equipment will tilt forward.

⚠ DANGER OF SQUEEZING FINGERS