



MOUNTING REFERENCE GUIDE

Monitoring and PDMS Mounting Solutions for GE Healthcare Carestation 750

CIM Cable Integrated Mounting System



THINK BEFORE
YOU PRINT



Content

General notes	3
Product highlights	4
How to create part numbers	5
Load capacity support arm and monitor adaptation	7
Adapter for GE Carestation 750	8
Monitor mount	
Extended pivot arm heavy load	9
Height adjustable arm	11
Extended height adjustable arm	13
Bracket and keyboard tray	
Extended pivot arm heavy load with W-bracket	15
Short height adjustable arm with W-bracket	17
Counter top	19

General notes

Product safety

CIM med® medical mounting solutions for anaesthesia machines are designed and tested to the most exacting of standards, are CE marked and adhere to the regulations (EU) 2017/745 (MDR). EN 60601-1 compliancy forms the basis of everything we do, to ensure that you and your customer get the safest, most beneficial solutions available in the world today.

Each arm is earthed back to the carts own earthing system. The arms can be restricted in movements to protect elements of the cart like rotameters, screens and cylinder manifolds. This also helps tip test conformity.

We have taken the time and effort to carry out all other testing such as static and dynamic load testing in combination with your anaesthesia cart. This is an essential part of medical compliance. We consider it our responsibility to ensure that the CIM configuration offered in this mounting reference guide, is tested, stable and safe. We will never supply you with a solution that is not 100% compliant. It is important that the maximum load carrying capacity for the anaesthesia machine as stated by the manufacturer, is not exceeded when using multiple configurations.

Please ask for a variety of TÜV-tested mounting solutions for anaesthesia applications. A declaration of conformity as well as the TÜV report can be provided on request: enquiry@cim-med.com

Mounting positions

Monitoring and PDMS can only be mounted on the right side of the Carestation 750.

Product marking

MD A stand-alone product that is already classified as a medical product class 1 before being integrated into an overall medical device. The product is CE-marked. A declaration of conformity is available.

~~MD~~ Accessories or spare parts that have no independent function and are therefore not subject to the MDR. These products are not CE-marked. A declaration of conformity is not required.

Warranty

All CIM products have a 5-year warranty to be free of defects in material and functionality from the date of delivery.

Product highlights

Non-height adjustable arm

Pivot and extended pivot arms offer **full cable integration** and easy cleaning for **infection control** in all critical care areas in the clinical environment. Our **5-years warranty** guarantees reliable, trouble-free and maintenance-free operation. The swivel function at the fixing point as well as the **tilt/ swivel function** of the monitor adaptation offers optimum positioning of your device. The extended pivot arm is a two-part arm with central pivot point between the two arm sections. All bottom-mounted monitors have the safety feature of a **20° tilt stop**.

Height adjustable arm

The height adjustable arm offers **full cable integration**, one-wipe clean surfaces and our **5-year warranty** on the entirety of the product. Vertically and horizontally adjustable, the arm offers optimal positioning of your device. A constant viewing angle when vertically adjusted is guaranteed by a parallel linkage. **Weight-optimized springs counter-balance** the monitor without making any adjustments. All bottom-mounted monitors have the safety feature of a **20° tilt stop**. Individual rotational stops may be integrated in any joint.

To provide the perfect combination of arm and monitor, the CIM articulated height adjustable arm has a simple quarter-turn lock to hold the arm safely in position when removing the monitor. Furthermore the arm is equipped with the unique safety feature of a slow-release gas spring. Should anybody try to remove the monitor before locking the arm in position, the arm will only rise at a slow and even tempo.

SAFETY

Monitor mount

How to create part numbers



1. Choose the arm version required:

H33013	=	extended pivot arm heavy duty
431	=	height adjustable arm
73113	=	extended height adjustable arm

2. Choose the monitor adaptations required:

A	=	5" plunge plate
G	=	VESA 75/100
PH	=	Philips table top mount
RR	=	infusion pump mount (only with H33013)

3. Mounting side:

4S-04R	=	GE Carestation 750, right side
--------	---	--------------------------------

4. Choose the colour requested for decorative parts (C):

gr	=	RAL 7024 (graphite grey)
bl	=	RAL 5013 (cobalt blue)

5. Colour of the aluminium die-cast parts:

1	=	RAL 9016 (traffic white)
---	---	--------------------------

6. Specify the spring depending on the weight of your monitor/ device - only with height adjustable arms (X):

60	=	2 - 5 kg / 4.4 - 11 lbs
95	=	4 - 9 kg / 8.8 - 19.8 lbs
110	=	4.5 - 12.5 kg / 9.9 - 27.6 lbs
130	=	7.5 - 13 kg / 16.5 - 28.7 lbs
160	=	12.5 - 18 kg / 27.6 - 39.7 lbs
175	=	13 - 16 kg / 28.7 - 35.3 lbs
220	=	16 - 18 kg / 35.3 - 39.7 lbs

Example

Grounding		01 + 02 Arm version + Monitor adaption		03 Mounting side		04 + 05 Colour		06 Spring tension
GE	-	431A	-	4S-04R	.	C1	-	X
GE	-	431A	-	4S-04R	.	bl1	-	95

Bracket and keyboard tray

How to create part numbers

1. Choose the arm version required:

H33013F	=	extended pivot arm heavy duty with W-bracket
73113F	=	extended height adjustable arm with W-bracket

2. Mounting side:

4S-04R	=	GE Carestation 750, right side
--------	---	--------------------------------

3. Choose the colour requested for decorative parts (C):

gr	=	RAL 7024 (graphite grey)
bl	=	RAL 5013 (cobalt blue)

4. Colour of the aluminium die-cast parts:

1	=	RAL 9016 (traffic white)
---	---	--------------------------

5. Specify the spring depending on the weight of your monitor/ device - only with height adjustable arms (X):

95	=	2.5 - 5.5 kg / 5.5 - 12.1 lbs
110	=	3.5 - 9.5 kg / 7.7 - 20.9 lbs
130	=	6 - 10 kg / 13.2 - 22 lbs
160	=	8 - 14.5 kg / 17.6 - 32.0 lbs
175	=	11 - 14 kg / 24.3 - 30.9 lbs

6. Specify the requested size of the keyboard tray (T):





30	=	300 x 210 mm / 11.8 x 8.3"
40	=	400 x 150 mm / 15.7 x 5.9"
50	=	500 x 200 mm / 19.7 x 7.9"
58	=	580 x 200 mm / 22.8 x 7.9"





Example

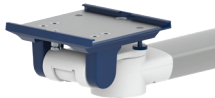Grounding		01 Arm version		02 Mounting side		03 + 04 Colour		05 Spring tension		06 Keyboard tray
GE	-	73113F	-	4S-04R	.	C1	-	X	-	T
GE	-	73113F	-	4S-04R	.	bl1	-	95	-	30

Load capacity support arm and monitor adaptation

	5" PLUNGE PLATE	VESA 75/100	PHILIPS TABLE TOP MOUNT	INFUSION PUMP MOUNT
				
RIGHT SIDE				
Arm version	H33013A	H33013G	H33013PH	H33013RR
Max. load capacity support arm	18.0 kg / 39.7 lbs	18.0 kg / 39.7 lbs	18.0 kg / 39.7 lbs	23.0 kg / 50.7 lbs
Max. load monitor adaptation	18.0 kg / 39.7 lbs	18.0 kg / 39.7 lbs	18.0 kg / 39.7 lbs	23.0 kg / 50.7 lbs
Product weight	3.3 kg / 7.3 lbs	3.4 kg / 7.5 lbs	3.2 kg / 7.1 lbs	2.8 kg / 6.2 lbs
Arm version	431A	431G	431PH	
Max. load capacity support arm	18.0 kg / 39.7 lbs	18.0 kg / 39.7 lbs	18.0 kg / 39.7 lbs	
Max. load monitor adaptation	18.0 kg / 39.7 lbs	18.0 kg / 39.7 lbs	18.0 kg / 39.7 lbs	
Product weight	4.4 kg / 9.8 lbs	4.3 kg / 9.5 lbs	4.3 kg / 9.5 lbs	
Arm version	73113A	73113G	73113PH	
Max. load capacity support arm	15.0 kg / 33.1 lbs	15.0 kg / 33.1 lbs	15.0 kg / 33.1 lbs	
Max. load monitor adaptation	15.0 kg / 33.1 lbs	15.0 kg / 33.1 lbs	15.0 kg / 33.1 lbs	
Product weight	5.2 kg / 11.5 lbs	5.1 kg / 11.3 lbs	5.0 kg / 11.0 lbs	

	W-BRACKET	W-BRACKET
		
RIGHT SIDE		
Arm version	H33013F	73113F
Max. load capacity support arm	18 kg / 39.7 lbs	14.5 kg / 32.0 lbs
Max. load monitor adaptation	18 kg / 39.7 lbs	14.5 kg / 32.0 lbs
Product weight *	5.0 - 5.5 kg / 11.0 - 12.1 lbs	6.7 - 7.1 kg / 14.8 - 15.7 lbs

*) Product weight varies depending on the size of the keyboard tray.

Keyboard with fold-up function available on request. The overall weight capacity of the fold-up version may be less due to the higher weight of the Z-bracket and fold-up keyboard tray.

Adapter for GE Carestation 750

Adapter for GE (DX-AES-01)

The adapter shown is mounted to the side channel of anaesthesia machine and has to be ordered separately.

Product weight: 0.5 kg / 1.1 lbs



Monitor mount

Extended pivot arm heavy duty

Monitor adaption, weights & loads

see page 7

Mounting side

GE Carestation 750, right side

Colours to choose **C**

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white

Part numbers

GE-H33013A-4S-04R.C1	5" plunge plate
GE-H33013G-4S-04R.C1	VESA 75/100
GE-H33013PH-4S-04R.C1	Philips table top mount
GE-H33013RR-4S-04R.1	infusion pump mount

NOTE: The arm can only be mounted on the right side of the anaesthesia machine.



MD




Illustration shown with 5" plunge plate

Monitor mount


Extended pivot arm heavy duty

Technical data

FRONT VIEW

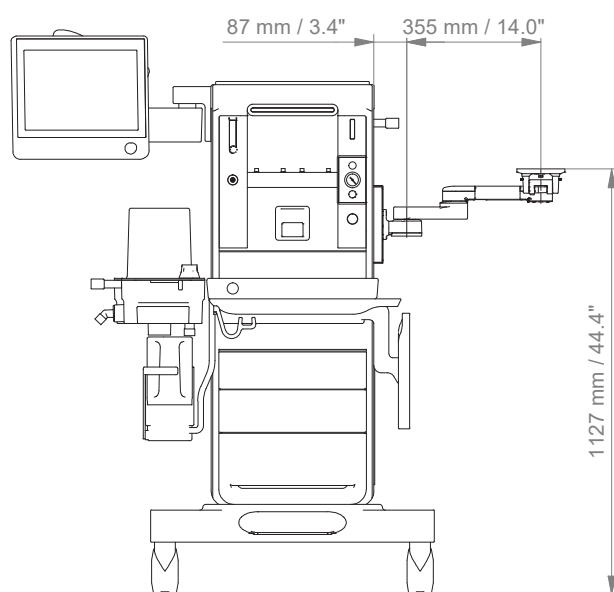

TOP VIEW

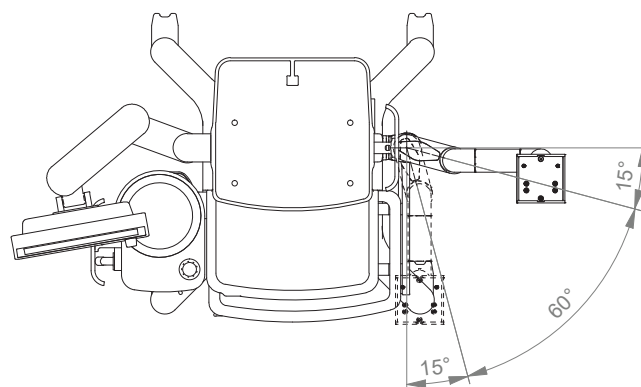

Monitor mount

Height adjustable arm

Monitor adaption, weights & loads

see page 7

Mounting side

GE Carestation 750, right side

Colours to choose **C**

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white

Height adjustment

slow-flowing spring, parallel linkage for a constant viewing angle

Balancing range **X**

600 N (2 - 5 kg / 4.4 - 11 lbs); 950 N (4 - 9 kg / 8.8 - 19.8 lbs); 1100 N (4.5 - 12.5 kg / 9.9 - 27.6 lbs);
1300 N (7.5 - 13 kg / 16.5 - 28.7 lbs); 1600 N (12.5 - 18 kg / 27.6 - 39.7 lbs);
1750 N (13 - 16 kg / 28.7 - 35.3 lbs); 2200 N (16 - 18 kg / 35.3 - 39.7 lbs)

Part numbers

GE-431A-4S-04R.C1_X	5" plunge plate
GE-431G-4S-04R.C1_X	VESA 75/100
GE-431PH-4S-04R.C1_X	Philips table top mount

NOTE: The arm can only be mounted on the right side of the anaesthesia machine.



MD




Illustration shown with 5" plunge plate

Monitor mount


Height adjustable arm

Technical data

FRONT VIEW

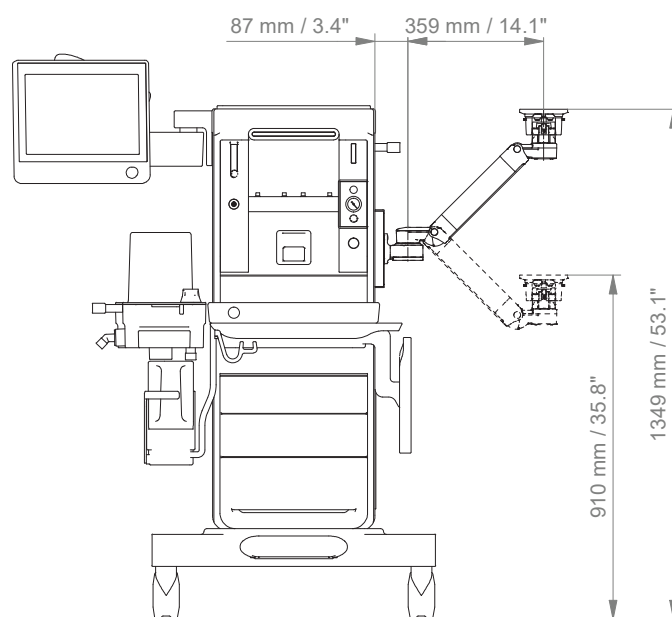

TOP VIEW

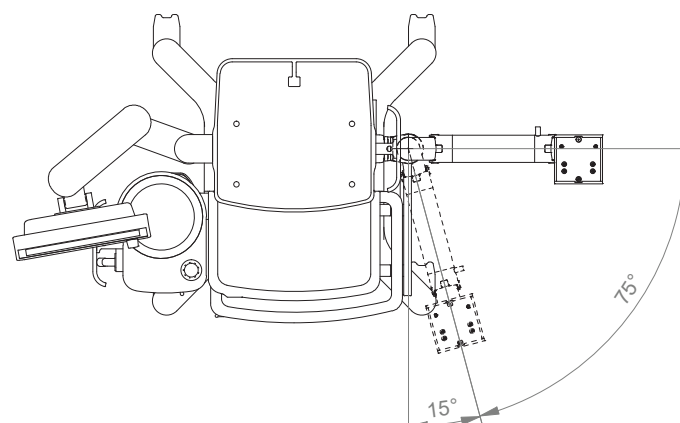

Monitor mount

Extended height adjustable arm



MD

Monitor adaption, weights & loads

see page 7

Mounting side

GE Carestation 750, right side

Colours to choose C

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white

Height adjustment

slow-flowing spring, parallel linkage for a constant viewing angle

Balancing range X

600 N (2 - 5 kg / 4.4 - 11 lbs); 950 N (4 - 9 kg / 8.8 - 19.8 lbs); 1100 N (4.5 - 12.5 kg / 9.9 - 27.6 lbs);
1300 N (7.5 - 13 kg / 16.5 - 28.7 lbs); 1600 N (12.5 - 15 kg / 27.6 - 33.1 lbs);
1750 N (13 - 15 kg / 28.7 - 33.1 lbs)

Part numbers

GE-73113A-4S-04R.C1_X 5" plunge plate
GE-73113G-4S-04R.C1_X VESA 75/100
GE-73113PH-4S-04R.C1_X Philips table top mount

NOTE: The arm can only be mounted on the right side of the anaesthesia machine.

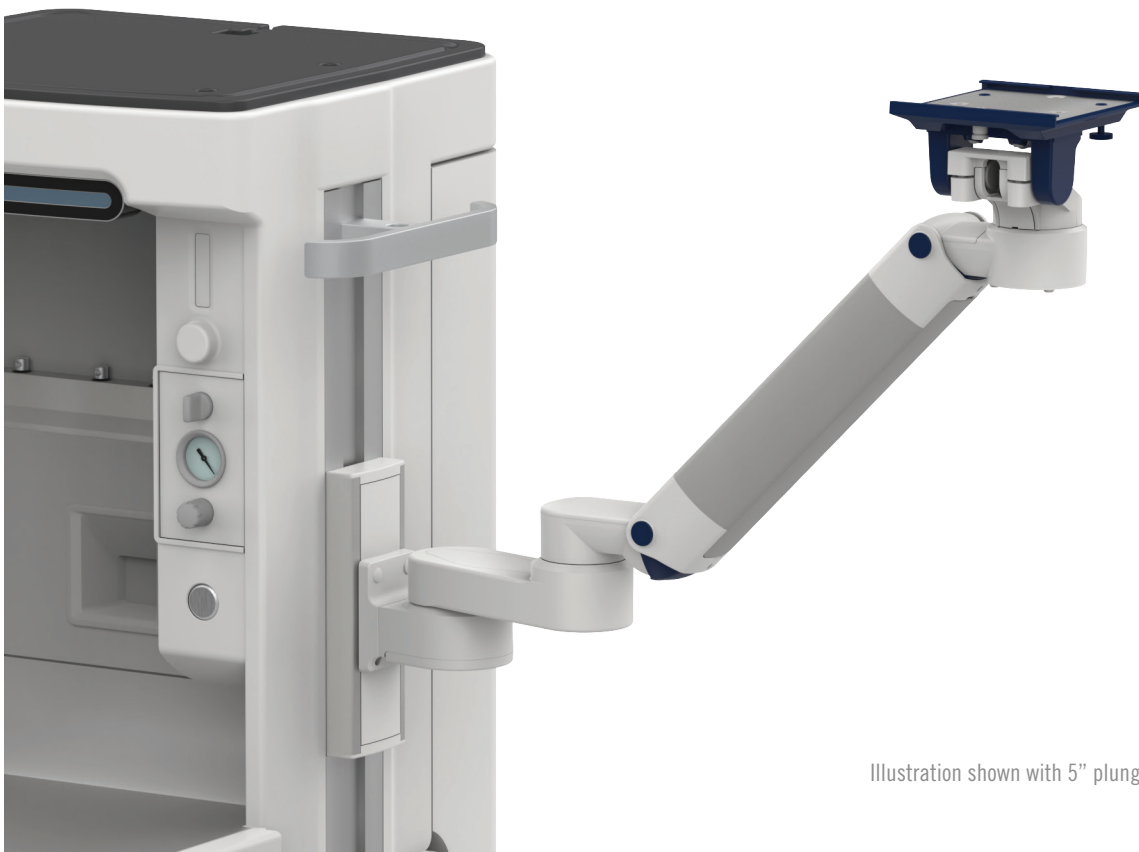



Illustration shown with 5" plunge plate

Monitor mount


Extended height adjustable arm

Technical data

FRONT VIEW

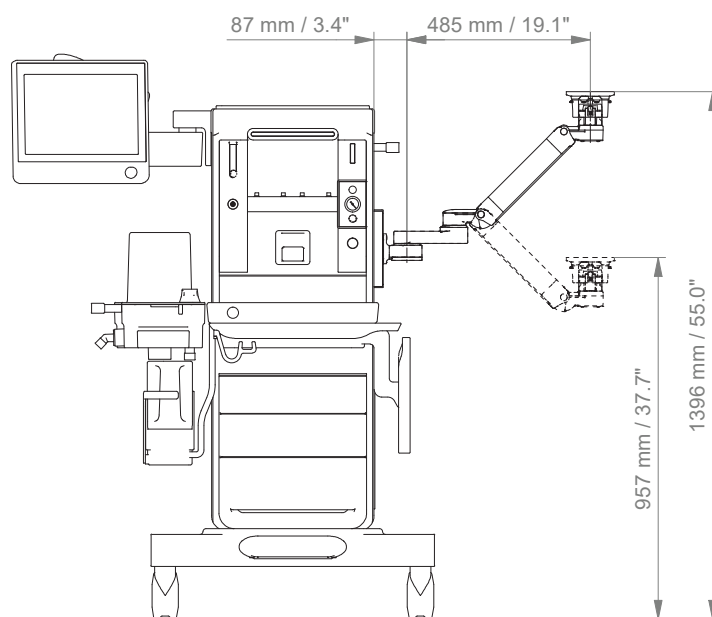

TOP VIEW

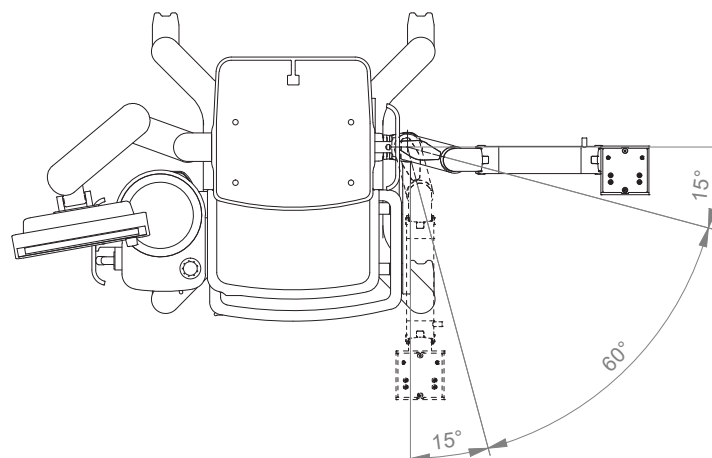

Bracket and keyboard tray

Extended pivot arm heavy duty with W-bracket

Monitor/ keyboard adaption

W-bracket VESA 75/100 and keyboard / mouse tray bracket with keyboard fold-up funktion on request

Weights & loads

see page 7

Mounting side

GE Carestation 750, right side

Colours to choose **C**

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white

Keyboard tray dimension **T**

300 x 210 mm, 400 x 150 mm, 500 x 200 mm, 580 x 200 mm
11.8 x 8.3", 15.7 x 5.9", 19.7 x 7.9", 22.8 x 7.9"


Part number

GE-H33013F-4S-04R.C1_T

NOTE: The arm can only be mounted on the right side of the anaesthesia machine.



MD

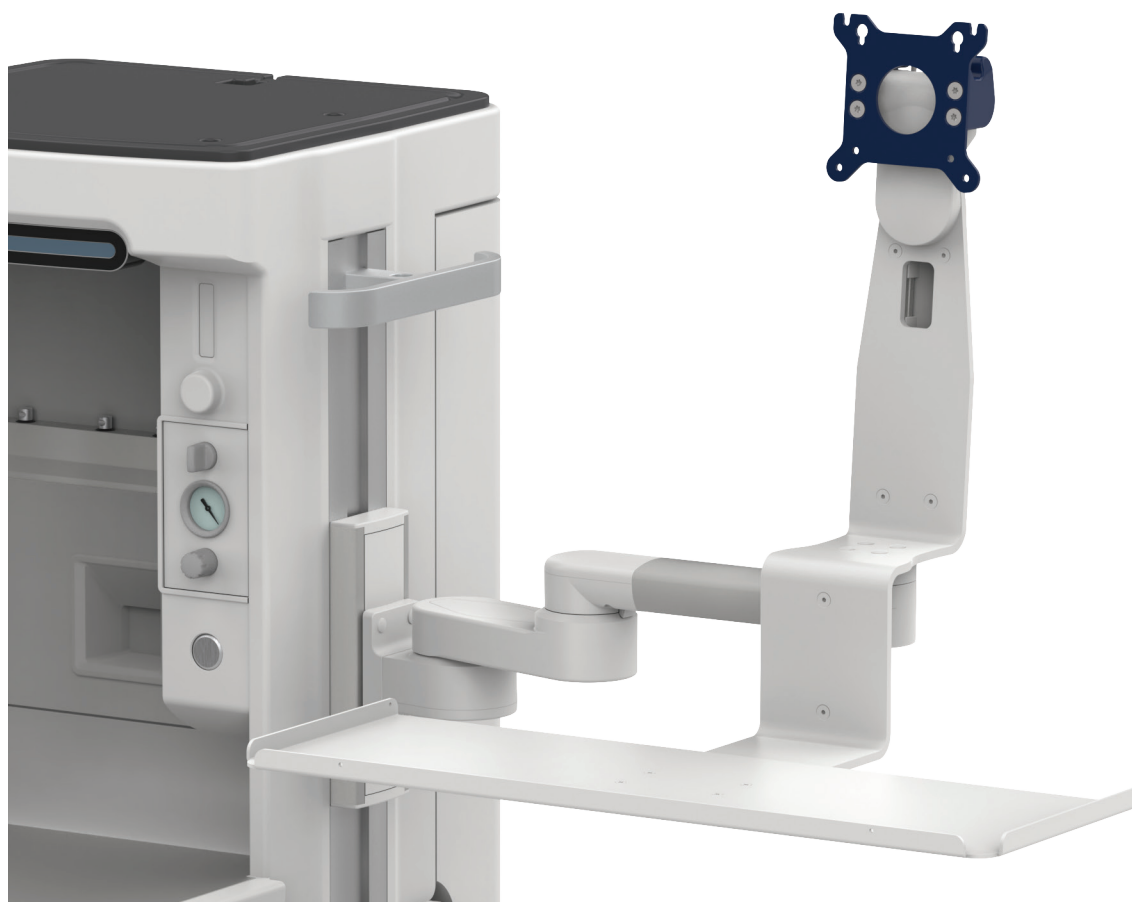


Bracket and keyboard tray


Extended pivot arm heavy duty with W-bracket

Technical data

FRONT VIEW

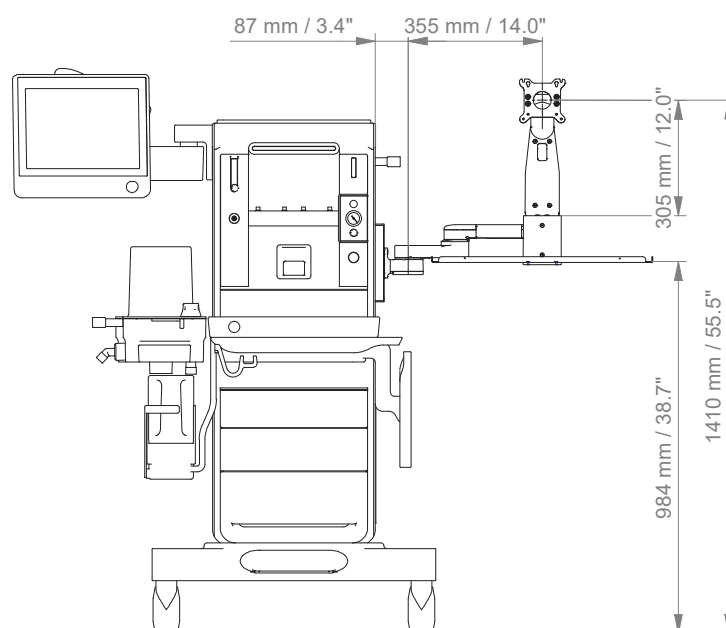

TOP VIEW

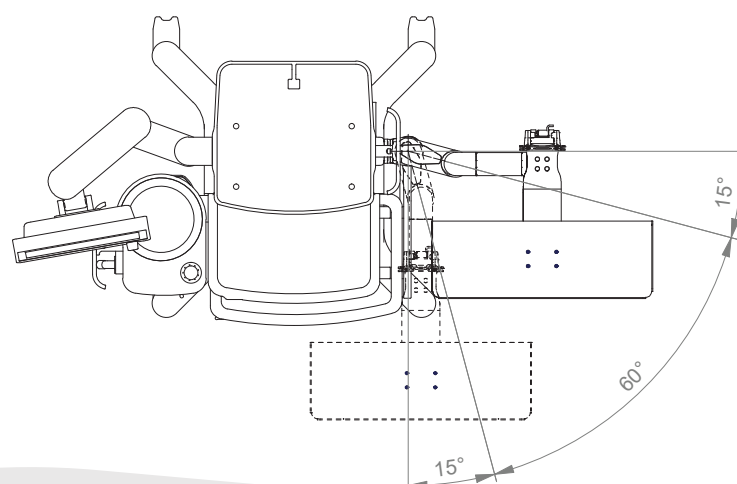

Bracket and keyboard tray

Extended height adjustable arm with W-bracket



Monitor/ keyboard adaption

W-bracket VESA 75/100 and keyboard / mouse tray bracket with keyboard fold-up funktion on request

Weights & loads

see page 7

Mounting side

GE Carestation 750, right side

Colours to choose **C**

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey
aluminum parts RAL 9016 traffic white

Balancing range **X**

950 N (2.5 - 5.5 kg / 5.5 - 12.1 lbs); 1100 N (3.5 - 9.5 kg / 7.7 - 20.9 lbs);
1300 N (6 - 10 kg / 13.2 - 22 lbs); 1600 N (8 - 14.5 kg / 17.6 - 32.0 lbs);
1750 N (11 - 14 kg / 24.3 - 30.9 lbs)

Keyboard tray dimension **T**

300 x 210 mm, 400 x 150 mm, 500 x 200 mm, 580 x 200 mm
11.8 x 8.3", 15.7 x 5.9", 19.7 x 7.9", 22.8 x 7.9"

Part number

GE-73113F-4S-04R.**C1_X_T**

NOTE: The arm can only be mounted on the right side of the anaesthesia machine.

MD




Bracket and keyboard tray


Extended height adjustable arm with W-bracket

Technical data

FRONT VIEW

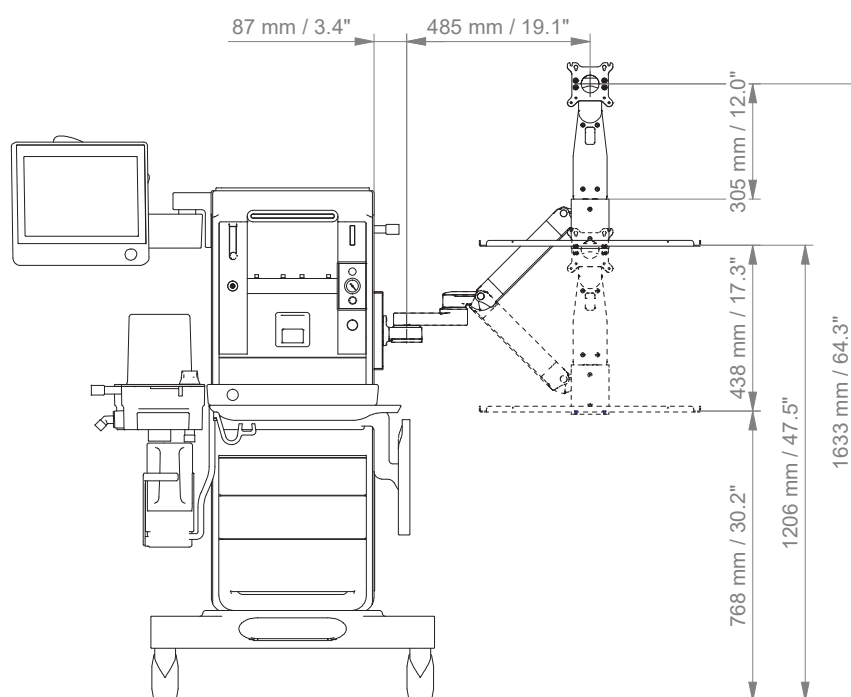

TOP VIEW

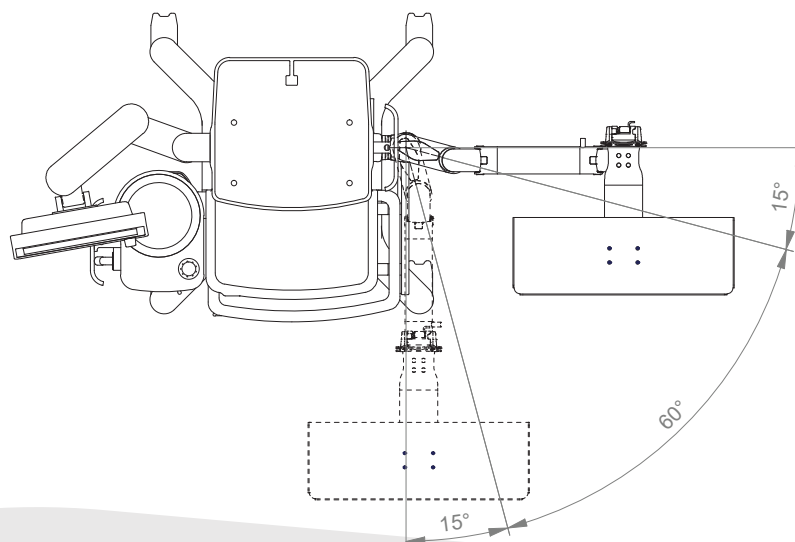

Top: Monitor mount

Counter top

Max. load monitor adaptation 18.0 kg / 39.7 lbs

Product weight 1.8 kg / 4.0 lbs

Colours to choose C

decorative parts RAL 5013 cobalt blue / RAL 7024 graphite grey

aluminum parts RAL 9016 traffic white

Part numbers

GE-HA-CS.C1_50

5" plunge plate

GE-HPH-CS.C1_50

Philips table top mount



MD

Illustration shown with 5" plunge plate




Top: Monitor mount


Counter top

Technical data

FRONT VIEW

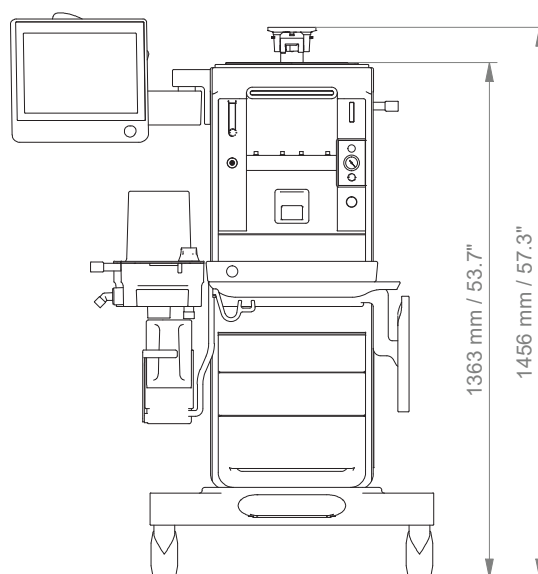

TOP VIEW

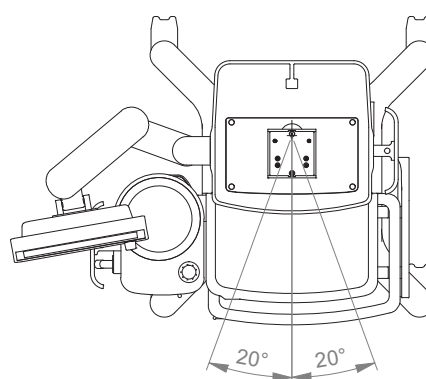



CIM med GmbH

Euro-Industriepark
Margot-Kalinke-Str. 9
80939 Munich / Germany

T +49 89. 978 94 08-0

CIM med Inc.

University Executive Park
301 McCullough Drive, Suite 400
Charlotte NC, 28262 / USA

T +1 704 221 0721

info@cim-med.com
www.cim-med.com

